

The logo for Dynamic Dental Equipment is centered in a red rectangular box. It features the word "DYNAMIC" in a bold, white, sans-serif font, with a stylized crosshair symbol integrated into the letter 'A'. Below "DYNAMIC" is the phrase "DENTAL EQUIPMENT" in a smaller, white, sans-serif font. The background of the entire image is a dark, green-tinted photograph of a dental procedure, showing a patient's head in a chair and a dental handpiece being used on a tooth.

DYNAMIC
DENTAL EQUIPMENT


REPLACING DYNAMIC DISC GUIDE

YOU WILL NEED:

AVAILABLE FROM YOUR DYNAMIC DISTRIBUTOR:

- DYNAMIC DISC REMOVAL TOOL
- DYNAMIC HEAD LOCKING TOOL
- NEW DYNAMIC REPLACEMENT DISC

- REMOVE HEAD FROM HAND
PIECE
- SCREW THE HEAD LOCKING
TOOL ON TO THE BACK OF THE
HEAD CASING

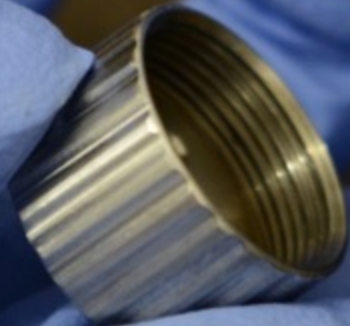

▪ PLACE THE DISC REMOVAL
TOOL OVER THE DISC

▪ ENGAGE THE
DISC REMOVAL
TOOL OVER THE
DISC AND
ROTATE ANTI-
CLOCKWISE

▪ THEN
REMOVE THE
OLD DISC AND
DISCARD



- INSERT THE NEW DISC INTO THE DISC REMOVAL TOOL AND PLACE INTO THE GUARD. TIGHTEN THE DISC IN A CLOCK-WISE ROTATION UNTIL SECURE. REMOVE LOCKING TOOL AND SECURE HEAD BACK ON TO THE HAND PIECE

